




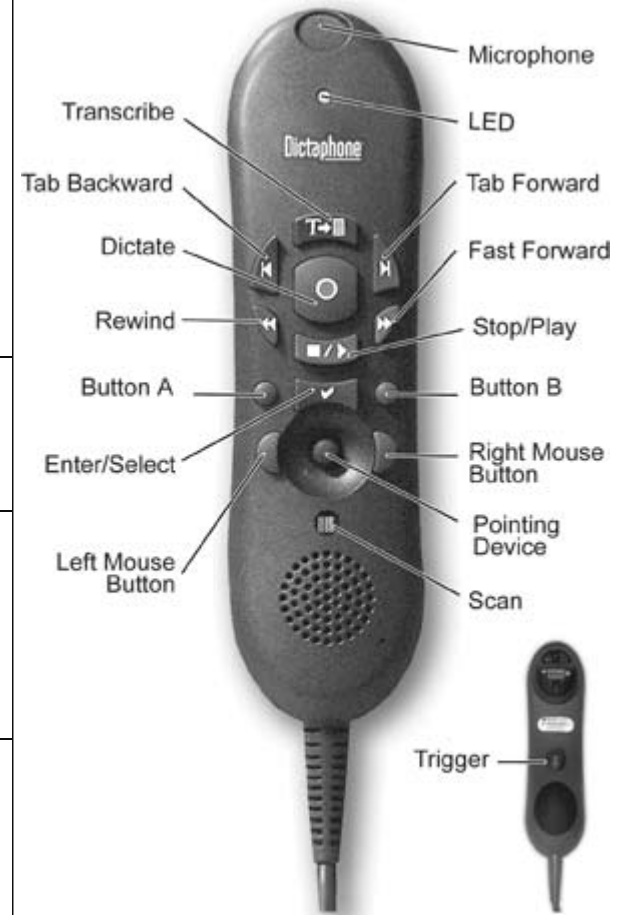









## Using the Dictaphone PowerMic II with Dragon NaturallySpeaking


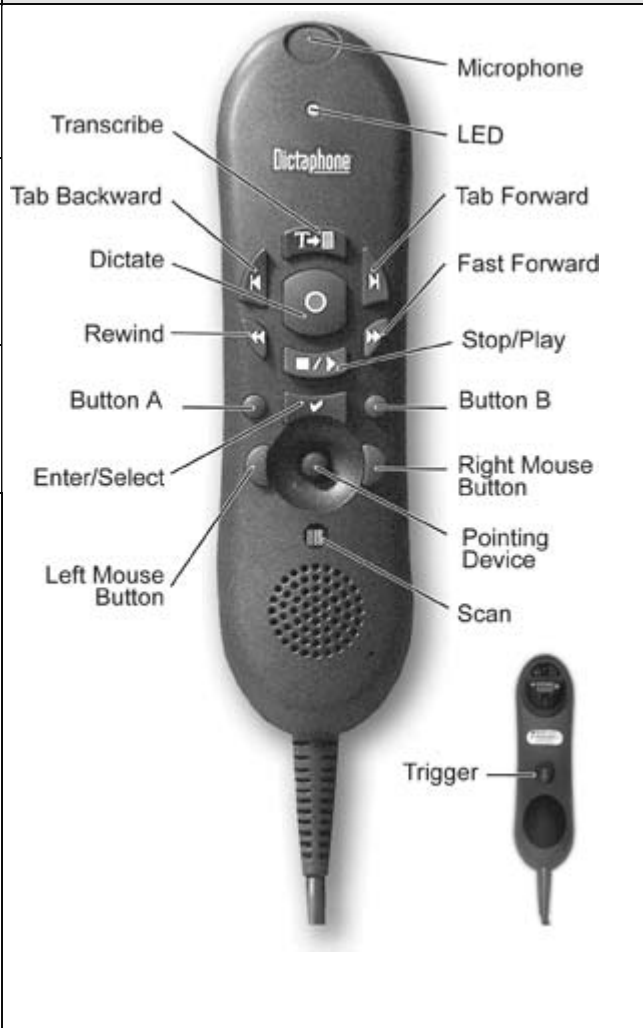



Dragon NaturallySpeaking Medical Edition has built-in support for the black Dictaphone PowerMic II.

The table below describes the black PowerMic II functionality with Dragon NaturallySpeaking and DragonConnect. For more information, see the online help.

Button	Name	In Dragon NaturallySpeaking	In DragonConnect
	LED	<p>The LED turns on when the user:</p> <ul style="list-style-type: none"> <li>Presses the Dictate button on the PowerMic</li> <li>Pressed the plus key (+) on the numeric keypad.</li> <li>Clicks the microphone button on the DragonBar.</li> <li>Clicks the microphone icon in the Windows taskbar.</li> </ul>	The same as defined in Dragon NaturallySpeaking.
	Transcribe	Press to end your current dictation and open the DragonConnect window when available.	Send file to the Enterprise Express.
	Dictate	<p>Turns the microphone on and off:</p> <ul style="list-style-type: none"> <li>Press to begin dictating your report.</li> <li>Press a second time to stop the dictating</li> </ul>	The same as defined in Dragon NaturallySpeaking.
	Tab Backward	Opens the Correction menu.	The same as defined in Dragon NaturallySpeaking.
	Tab Forward	None	Press this key to execute the keyboard tab command to move from field to field.



Button	Name	In Dragon NaturallySpeaking	In DragonConnect	
	Rewind	None	Press this key to: <ul style="list-style-type: none"> <li>• Scroll up a drop-down list.</li> <li>• Rewind word by word in an edit field.</li> </ul>	
	Fast Forward	None	Press this key to: <ul style="list-style-type: none"> <li>• Scrolls down a drop-down list.</li> <li>• Fast Forward word by word in an edit field.</li> </ul>	
	Stop/Play	Playback This hotkey causes Dragon NaturallySpeaking to play back your dictation so you can more easily identify what you meant to say during correction.	The same as defined in Dragon NaturallySpeaking.	
	Enter/Select	None	Press to activate the keyboard enter command.	
	Button A: Left Customizable button	Force command recognition. This hot key forces Dragon NaturallySpeaking to recognize what you say as a command..	The same as defined in Dragon NaturallySpeaking.	
	Button B: Right Customizable button	Force dictation recognition. This hot key forces Dragon NaturallySpeaking to recognize what you say as dictation.	The same as defined in Dragon NaturallySpeaking.	

Button	Name	In Dragon NaturallySpeaking	In DragonConnect	
	Left mouse/right mouse button	Functions in the same way of the buttons of the mouse	The same as defined in Dragon NaturallySpeaking.	
	Pointing device	Used to position the mouse insertion point and select objects as would be done with the mouse	The same as defined in Dragon NaturallySpeaking.	
	Scan	For microphones with a scanner: activates scanner.	The same as defined in Dragon NaturallySpeaking.	
	Trigger (on rear of microphone)	Same function as left mouse button.	The same as defined in Dragon NaturallySpeaking.	